

UNIVERSITY OF PITTSBURGH

# STEPUP Observing Guidelines



*Allegheny Observatory*

*Author:* STEPUP Team

*Supervisor:*  
Professor Wood-Vasey

# Contents

<b>1</b>	<b>Opening Checklist</b>	<b>2</b>
1.1	Checking the Conditions . . . . .	2
1.2	Preliminary Setup . . . . .	2
1.3	General Startup . . . . .	3
1.4	Calibration . . . . .	3
1.5	Focusing . . . . .	5
1.6	Guiding . . . . .	6
1.7	Sequencing . . . . .	6
<b>2</b>	<b>Abridged Opening Checklist</b>	<b>7</b>
<b>3</b>	<b>Closing Checklist</b>	<b>9</b>
<b>4</b>	<b>Troubleshooting</b>	<b>10</b>
4.1	WinTV . . . . .	10
4.2	No stars in Image . . . . .	10
4.3	Wrong star field . . . . .	10
4.4	Dome . . . . .	11
4.5	Not focusing . . . . .	11
<b>5</b>	<b>Appendix</b>	<b>12</b>
5.1	Sample Observing Report . . . . .	12

# Chapter 1

## Opening Checklist

Chapters 1 and 2 provide instructions on the complete observing process for STEPUP purposes. The following instructions do not require any knowledge of physics or astronomy. Therefore, anyone who can read and follow instructions can observe. If you are an experienced observer, please go to Chapter 2 for an abridged opening checklist. Please refer to Chapter 4 for troubleshooting information.

### 1.1 Checking the Conditions

- Use Pittsburgh Clear Sky Chart (<http://cleardarksky.com/c/PittsburghPAkey.html>) during the days before to check the possibility of observing that night.
- The day of, check NOAA's IR maps (<http://www.goes.noaa.gov/ECIR4.html>) to look at cloud cover in addition to Pittsburgh Clear Dark Sky.
- At 7:15pm you can start checking the All Sky Camera on Allegheny Observatory's website (<http://www.pitt.edu/~aobsvtry/All-Sky.html>) to directly verify the cloud coverage at the observatory. Looking out the window yourself has also shown to be a good way of checking for clouds.

### 1.2 Preliminary Setup

- In Allen 302, log on to the Ptolemy computer and open up the terminal.
- Create a directory for the target you are observing in the raw folder (`/home/depot/STEPUP/raw`) if one does not exist already. Within this newly created directory, make a directory for the data you will take that day in the format of year-month-day. For example, if you are observing the target HAT-P-13 on November 12, 2013 and a directory already exists for that target within raw, you should type `mkdir /home/depot/STEPUP/raw/HAT-P-13/2013-11-12`.
- Copy the template for the observing report into the directory that you just created by typing the command `cp /home/depot/STEPUP/obsreport-template.txt` (type a space) `/home/depot/STEPUP/raw/target/date`. You will type observing report here by editing this template. See the appendix for a sample observing report.
- On the Lamashtu computer, open the terminal and type: `vncviewer 136.142.17.99`. A prompt will ask you for a password, which is **foo**. You want to record four pieces of information

in the observing report: *cloud coverage, temperature, wind, and humidity*. They are located on the ClarityII window while the seeing camera is shown as a plot.

GitHub pushing test

### 1.3 General Startup


1. On the Ptolemy computer, open the terminal and type: **vncviewer ao-keeler.phyast.pitt.edu**. A prompt will ask you for a password, which is **foo**. If the computer is not logged on, you will have to log on with the password: **group1**. The username should be defaulted as *group1*, but if it isn't, change the username to *group1*. The username and password are the same.
2. Open the **Logitech QuickCam** on the desktop. When a window pops up, click on the first icon in that window. Use this to monitor the movement of the telescope and dome.
3. Open the **Keeler Power Controller 1.1** on the desktop. You will see a series of on and off buttons. Turn the telescope on first (an automated one minute waiting period follows) followed by the rest of the controls. After you have turned on the **Lens Cap Power**, click on the *Open* button under **Lens Cap**. View the webcam to make sure the lens cap is opening.
4. Open the **Digital Dome Works** on the desktop. Two windows should pop up: *Digital DomeWorks v5.2* and *RCX Control*. **WARNING: KEEP THE RCX WINDOW ON LEFT SIDE OF SCREEN. IT MUST NOT MOVE TO THE RIGHT HALF OF THE SCREEN!** If the RCX window even briefly edges over to right side of screen, call Lou.
5. Open the **MaxIM DL** on the desktop. When the MaxIM window pops up, click on the **Toggle CCD Control** button in the row near the top of the window. When the MaxIM CCD window opens, you should be in the *Setup* tab. Click on the *Connect* button followed by the *Cooler On* button.
6. Still in the MaxIM DL window, click on the **Toggle Telescope Control** button in the same row near the top of the window. When the Telescope Control window opens, you should be in the *Setup* tab. Click *Connect* for both the telescope and focuser.

### 1.4 Calibration

You should be taking calibrations each night that you observe. However, if you have time constraints (For example, if a transit starts at 9pm, but you just started setting up at 8:15pm.), you can skip this section and save it for the end of the night.

1. In the Digital DomeWorks window, click on *Go To* and enter 270 when it prompts you to. The dome will likely stop shy of 270 degrees, so you will need to over shoot and click *STOP* a few degrees before it reaches 270. You may have to play with it a few times before the dome is positioned correctly. Use the webcam to verify the dome is moving.
2. In the RCX Control window, turn tracking off by clicking on the green button that says *Tracking* near the top. If there is a red button labeled *Not Tracking*, then it is already set. Regardless, the final setting should be *Not Tracking*.

3. In the RCX Control window, note the *ALT* and *AZ* values in the top right corner of the window. You want to set these to **ALT: 10 degrees** and **AZ: 180 degrees**. Within one degree is close enough. To do this, you will use the buttons in the *Motion Control* box in the center of the RCX Control window. Use 'N' and 'S' buttons to adjust the ALT value larger and smaller respectively. Use the 'W' and 'E' to adjust the AZ value larger and smaller respectively. The radio buttons below labeled *Guide*, *Center*, *Find*, *Slew* represent the magnitude of the response with *Slew* being the greatest and *Guide* the smallest. For example, if the ALT is at 40 deg, you will want to use *Slew*, but if it is at 11 deg, you will want to use *Guide*. If the direction arrows do not work, select *Guide* first and try. This will wake the program up if frozen.
4. You will now create a folder to store the nights images. Click on the Start menu -> My Computer -> group 1 on 'aoserver1\data'(Z:) -> STEPUPDataFiles. In this folder, create a new folder in the format of YYYY-MM-DD for the night of observation.
5. In the MaxIM CCD sub-window within the MaxIM DL window, click on the *Sequence* tab. Click on *Options* and select *Set Destination Path...* Select the folder you just created.
6. Click on *Options* again and select the *Setup Sequence* choice. Enable the two boxes Flat and Bias, setting the suffix as F and B respectively. Make sure the *Binning* is set to 1 and *Repeat* is set to 11 for each. No other settings should be enabled at this point. **IMPORTANT:** The bias images are instantaneous. This should be on default so the *Exposure* cannot be changed. The flat images have exposure times of 2 seconds. Biases will not have a filter, however, you should be sure to set the filter for the flat sequence to be the same as the filter you plan to take target images with. Indicate these filters in the suffix box. Your setup should look like the screenshot below. If you plan to take target images in several filters, you will need to, for each filter, set up a separate flat sequence with the same specifications as above in each different filter. Typically, for observing exoplanet transits, you can use the r filter. For variable stars, use B and V filters. For example, if you were observing a variable star your *Setup Sequence* window will look like this:

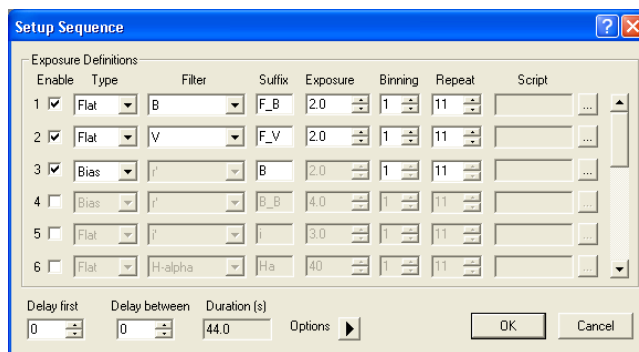

7. Still in the *Sequence* tab, enter the name of the star in the *Autosave Filename* box. Make sure that "1" is selected in the *Start at* box.
8. Start taking the calibration images by clicking *Start*. As the images begin to download, open the folder you created. Make sure the calibration images are saving to that folder. You also want to make sure that the flat field lights are on.

Side-note: We also take 30 minute dark calibration images, but this only needs to be done once a month and is usually not done on an observing night due to the fact that it takes over 1.5 hours. To take the 30 min darks the telescope and dome can be left where they are. In the Sequence tab instead of selecting Flat and Bias select Dark and put the exposure time at 1800 seconds (30 min). Other than that the procedure is the same.

## 1.5 Focusing

1. In Digital DomeWorks, move the dome back to 180 degrees using the *Go To* setting. When it nears 180 degrees, click *Home*. It should say that "Dome is Home" and "DSR closed". If not, continue to click *Home* until it says "Dome is Home"
2. Still in Digital DomeWorks, click the *Open* tab. This will open the dome shutters. Use the webcam to visually verify that the shutters of the dome are in fact opening.
3. Open **Starry Night College 6** on the desktop (If it asks for updates, do not update. If it asks for downloads, do not download). When the window opens click on the *Telescope* tab on the left side of the screen. Then click *Connect*. Minimize the *Telescope* tab. Do not use these tabs again until closing.
4. In Digital DomeWorks, check the *Slave to Telescope* box in the lower left corner.
5. Open **STVRemote** on the desktop. Click on the *Link* tab at the top and select *Establish Link (Com 3)*. Click on the *Image* button once. Using the buttons under the *Value* button, adjust the exposure time to 1.5 seconds. Finally, click the *Monitor* button twice.
6. In Starry Night College 6, click on the *Options* tab on the top menu bar. Go to *Orientation* and select *Equatorial*. Go to the *Edit* tab on the top menu bar and select *Centre On...* Enter the target coordinates.
7. Zoom in to the small rectangular box that should be labeled "Meade RCX-400..." The zoom controls are in the upper right corner of the window. Find a 6-8 magnitude star nearby that is not a binary star system by moving your cursor over the stars. Right click on the star and select *slew telescope here*. You should see a red cross hair cursor move to that point. This cursor indicates where the telescope is pointing.
8. Visually check on the webcam monitor that the telescope is pointing through the opening in the dome. Turn off flat field lights in the Keeler Power Controller 1.1.
9. In MaxIM CCD window, go to the *Focus* tab. Click on the *CCD* radio button. **Make sure binning is set to 3 when CCD radio button is selected!** Set the exposure time between 1-5 seconds depending on quality of image. *Delay* should be set at 0. Click *Start Focus*. A window of the field of view will open. Click *Stop*. The dimensions of this window should be the same as the dimensions of the CCD window: 3072 x 2048. Find the focus star in that image, right click on it and select *point telescope here*.
10. Go back to Starry Night College 6. The red cursor may or may not have shifted a bit. Right click on the focus star and select *sync on 'starname'*.
11. Open **FocusMax**, which should be in the taskbar (It opened automatically when you connected the focuser). Click *Focus*. Record the **HFD** that FocusMax returns in the observing

report. Also, under the *Inspect* tab of MaxIM CCD, record both the **FWHM** and **HFD** values.

## 1.6 Guiding

1. In Starry Night College 6, re-enter the target coordinates by going to *Edit -> Centre On...* Zoom into Meade RCX-400 box. Find a 6-8 magnitude star in the field to be your guide star. This star does not need to be the same as the focus star, and the focus star does not need to be in the Meade RCX-400 box when guiding. Move the box around until the guide star is located in the guide star box (the smaller box above) **AND** the target star is still within the Meade RCX-400 box. Do not rotate the guide star box. Ideally, we do not want the target star too far from the center. Right click on the center of the Meade RCX-400 box (not the target!) and select *slew telescope here*.
2. In the MaxIM CCD *Focus* tab, click the *Guider* radio button and set the exposure time to 3 seconds. Click *Start Focus*. The guide star should show up in the image that pops up. Click *Stop*. *\*\*Note: If the guide star does not show up, see Chapter 4: Troubleshooting*
3. Still in MaxIM CCD, go to the *Guide* tab and select the *Expose* radio button. Make sure the exposure time is set to 3 seconds and binning is set to 3. Click *Start*. An image with the guide star should pop up. Click on the guide star in the image. Back in the *Guide* tab, select the *Track* radio button and click *Start*. A very small image window with the guide star should show up. Monitor this window occasionally to ensure the guide star is still present and appears in focus.

## 1.7 Sequencing

1. In the MaxIM CCD *Sequence* tab, click on *Options -> Setup Sequence*. Uncheck the calibration images. Check a box where *Light* is selected. Select the appropriate filter and suffix (r' for exoplanets). Select an exposure time depending on the magnitude of the target star, normally between 30-120 seconds. Click *Ok*.
2. Back in the *Sequence* tab, make sure the image starts at 1 and the target star name is in the *Autosave Filename* box. Click *Start*.
3. After the image downloads, click on the cross hair button near the top of the MaxIM DL screen. A circular aperture now becomes your cursor, and a small window with information about the pixel counts should pop up. Move the circular cursor over the target star. If the counts are between 20,000-50,000, then you are set. Record the **FWHM**, **Maximum**, **Minimum**, **Bgd Avg** in the observing report. You will record this information every 30 minutes. If the counts are below 20,000, you will want to increase the exposure time. If the counts are above 50,000, you will want to decrease the exposure time.
4. If the maximum hits 65,535 this means the ccd has been maxed out. If the target star hits 65,535 in an image than we *cannot* use that image in our data analysis.
5. Try to keep the exposure time below 200s so we get a decent number of images of the the transit. If you observing a star dimmer than magnitude 11.5, this may mean that the maximum counts dip to around 10,000. That's okay, as long at its remains above 7,500 the data should be analyzable, just slightly noisier.

## Chapter 2

# Abridged Opening Checklist

1. Open all controls in the **Keeler Power Controller 1.1**.
2. Open Logitech webcam.
3. Open **Digital Dome Works**. Remember to keep **RCX** window on left side of the screen!!!
4. Open **MaxIM DL**. Open **Toggle CCD Control** and click both *Connect* and *Cooler On*. Open **Toggle Telescope Control** and connect both focuser and telescope.
5. Open **STVRemote**. Under the *Link* tab, click on *Establish Link (Com 3)*. Click *Image* button once. Set exposure time to 1.5 seconds. Click *Monitor* twice.
6. Create folder for the night in STEPUPDataFiles under the Z drive.
7. Set destination path to that folder in MaxIM CCD.
8. **Calibrations**
  - (a) Move telescope to **ALT: 10 deg** and **AZ: 180 deg** in the RCX window.
  - (b) Move dome to 270 deg.
  - (c) Ensure flat field lights are on.
  - (d) In MaxIM, setup the sequence with 11 images each of the flats and biases. Flats should be 2 seconds and biases are instantaneous. The suffixes are F and B.
  - (e) Make sure image starts at 1 and the filename is the name of the star you are observing.
  - (f) Click *Start*.
9. Move dome back to 180. Click *Home* and make sure DDW says it is home. Open the shutters. Click *Slave to Telescope*.
10. Open **Starry Night College 6**.
  - (a) Connect the telescope.
  - (b) Set orientation to *Equatorial*.
  - (c) Center onto the target coordinates. Right click on the field of view and select *Point Telescope Here*.



- (d) Select a star between the magnitudes 6-8 nearby.
11. Go to MaxIM DL. In the *Focus* tab, select CCD radio button and start 1-5 second focus. **Make sure binning is set to 3 when CCD radio button is selected!** Right click on target star and select *Point telescope here*.
  12. Back in Starry Night, right click on focus star and select *sync to 'starname'*.
  13. **Turn off flat field lights!**
  14. Go to FocusMax and click *Focus*. Record the **HFD** and **FWHM** info from both FocusMax and MaxIM.
  15. Re-enter coordinates in Starry Night. Put a guide star in the guide star box. Right click on center of field of view and select *Point telescope here*.
  16. In MaxIM Focus tab, select the *Guide* radio button and click *Start Focus*. Make sure the guide star is there.
  17. In MaxIM Guide tab, select the *Expose* radio button and click *Start*. Click on the guide star in the image that pops up. Back in the Guide tab, select the *Track* radio button and click *Start*.
  18. In MaxIM Sequence tab, setup the sequence so the light images are being taken. Make sure image starts at 1 and the *Autosave Filename* is filled in with target star name.

## Chapter 3

# Closing Checklist

1. Under the *Setup* tab in the MaxIM CCD window, press the *Warm Up* button. Wait about 1 minute. Click on *Cooler Off*. Click on *Disconnect*. Close the window.
2. In the MaxIM Telescope Control window, click *Disconnect* for both the telescope and focuser. Close the window.
3. Exit MaxIM DL completely.
4. In the Starry Night College 6 program, open the *Telescope* tab on the left. Disconnect the telescope. Exit Starry Night College 6, and don't click Save.
5. In the RCX Control window, click on the *Park* button. Don't close the window. When you close Digital Dome Works, the RCX Control window will close.
6. In Digital DomeWorks, unslave the telescope and turn the dome to 180 degrees by clicking GoTo. If the dome happens to get stuck, try turning it in the opposite direction. When the dome stops, click *Home* and make sure the message "Dome is Home" is there. Click *Close* and make sure the message "DSR is CLOSED" is there. Visually check the webcam screen to ensure dome is in fact closed (If it's dark, turn on flat field lights).
7. In STVRemote, click *Interrupt* and close the program.
8. In the Keeler Power Controller 1.1, click *Close* under the lens cap. Make sure lens cap closes in the webcam view.
9. Close the Logitech webcam.
10. Close Digital Dome Works.
11. Turn off all the power switches in the Keeler Power Controller 1.1 starting from *Lens Cap Power* up to *Telescope*. Close the window.
12. Log off the computer.

## Chapter 4

# Troubleshooting

### 4.1 WinTV

WinTV shows the image that the finder camera sees. The image is not as clear, but you should be able to see a few bright stars in WinTV. Match these stars with those in the images that you take. If the stars do not match up, the telescope is not pointed at where Starry Night says it is pointed.

- Click *start* -> *Remote Desktop Connection*. Log onto 136.142.17.105. The password is **blah**

### 4.2 No stars in Image

- Did you turn off the flat field lights? Use webcam to check.
- Are the shutters of the dome opened? Use webcam to check. If DDW says it is open, but the webcam shows it is not, move the dome back to 180. Click *Home*. Click *Close*. Click *Reset*. Click *Open*. If the dome still does not open, call Lou.
- What does the All Sky Camera show? Are there thick clouds?
- Check the sequence setup. Make sure we are taking images in R or r'.
- Did the dome get stuck? If so, *See Dome*.
- Is the dome still slaved to the telescope? If not, recheck that box in DDW.

### 4.3 Wrong star field

- *See WinTV*
- Did daylight savings recently happen? If so, call Lou.
- Go to Edit -> Centre On. Re-enter the target coordinates. Make sure the coordinates are correct before centering.

## 4.4 Dome

- If DDW shows the shutters are opened, but it is not, close the shutters in DDW. Click *reset*. Reopen the shutters. If it does not work, call Lou.
- If the dome is stuck, unslave the telescope. Move the dome in the opposite direction it was going by at least 90 degrees. Reslave the telescope. If the dome still gets stuck, move the dome in the opposite direction by at least 181 degrees. Then reslave the telescope. The dome should now reach where the telescope is from the other direction. *Note: In the winter, the dome sometimes gets stuck in the 300-360 degree range.*

## 4.5 Not focusing

- If you are having trouble focusing, try focusing on a magnitude 6ish star that is directly overhead (sidereal time equals star's RA) and a declination of 40 degrees.
- If you are unsure if the telescope is working (i.e. it is so out of focus you can't see anything) try pointing it at a really bright star such as Vega, Capella, Polaris or even Sirius. While you might be able to see these stars as bright blobs taking up half of the viewing window you won't be able to focus on them since they are so bright. However, they do serve as a good indicator that the telescope is working and pointing where you tell it to go.

# Chapter 5

## Appendix

### 5.1 Sample Observing Report

Observers: Anja, Helena, Meghan, Cynthia

Date: 2017-07-25

Scientific Goal of Observation: Contribute data for AAVSO Alert Notice 555 on WASP-65.

#### Conditions

Cloud Cover: Clear

Wind: 4.1mph

Humidity: 38%

Temperature: 35.6°F

The information below can be found by googling the target. Make sure the source is reliable. It is a good idea to cross reference sources to make sure the information is consistent.

Target: WASP-65

Mag: 11.9

RA: 08h 53m 17.38s

Dec: 08° 31' 22.8"

4:45pm - Began setup.

In this space below, record times at which you:

- Checked that dome and telescope were functioning.
- Began taking dark calibrations (if applicable).
- Began taking calibration images (flats and biases).
- Began focusing the telescope.

5:00pm - Visually confirmed that the dome and telescope are working.

5:05pm - Started taking dark calibrations.

Darks: 3                      Exposure Time: 1800s

7:30pm - Started taking flat and bias calibrations.

Flats: 11                      Filter(s): r  
Biases: 11

8:13pm - Focused telescope.

When you focus, FocusMAX will stop automatically when it focuses well. Record the values displayed in FocusMax and MaxIM here:

FocusMax Info -                      HFD: 4.79  
MaxIM Info -                      FWHM: 2.91                      HFD: 4.63

8:35pm - Began guiding.

Guide Star: TYC811-2091-1      Mag: 10.12  
RA: 8:53:50.3                      Dec: 8:43:01

8:40pm - Began sequencing.

Filter(s): r                      Exposure Time: 120s  
FWHM: 2.975                      Bkgr(avg): 4543  
Star(max): 20897                      Star(min): 4998

In this space, record image information every 30 min until you're finished observing with the following format:

While images are being taken, it is possible to lose focus, and refocusing will be necessary. Make notes between image checks of anything that needed to be redone or any environmental changes.

8:47pm -                      003r  
FWHM: 3.957                      Bkgr(avg): 4688  
Star(max): 21438                      Star(min): 4565

9:17pm -                      018r  
FWHM: 4.013                      Bkgr(avg): 5797  
Star(max): 23428                      Star(min): 4786

For image 19 we lost the guider and the image was completely blank. It looked like a dark calibration image. As of 9:10pm it seems to have recovered and all of our images look normal.

9:47pm -                      033r  
FWHM: 8.138                      Bkgr(avg): 5467  
Star(max): 18964                      Star(min): 3425

10:00pm - Refocused, but didn't seem to get any better. Images look okay though.

10:01pm - 039r

FWHM: 9.716 Bkgr(avg): 5347  
Star(max): 12134 Star(min): 2426

10:31pm - 054r

FWHM: 5.444 Bkgr(avg): 5023  
Star(max): 22452 Star(min): 2003

11:01pm - 069r

FWHM: 4.324 Bkgr(avg): 4897  
Star(max): 23982 Star(min): 3012

11:31pm - 084r

FWHM: 6.788 Bkgr(avg): 4563  
Star(max): 32094 Star(min): 2121

12:01am - 099r

FWHM: 7.323 Bkgr(avg): 5432  
Star(max): 23453 Star(min): 1890

12:31am - 114r

FWHM: 7.021 Bkgr(avg): 5567  
Star(max): 24221 Star(min): 2189

12:32am - Stopped sequence.

Science Images: 114 Filter(s): r

12:33am - Begin shutdown.

12:45am - Shut down.

Success: Yes